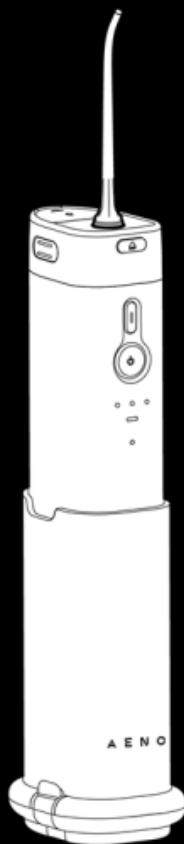
AENO

Quick Start Guide

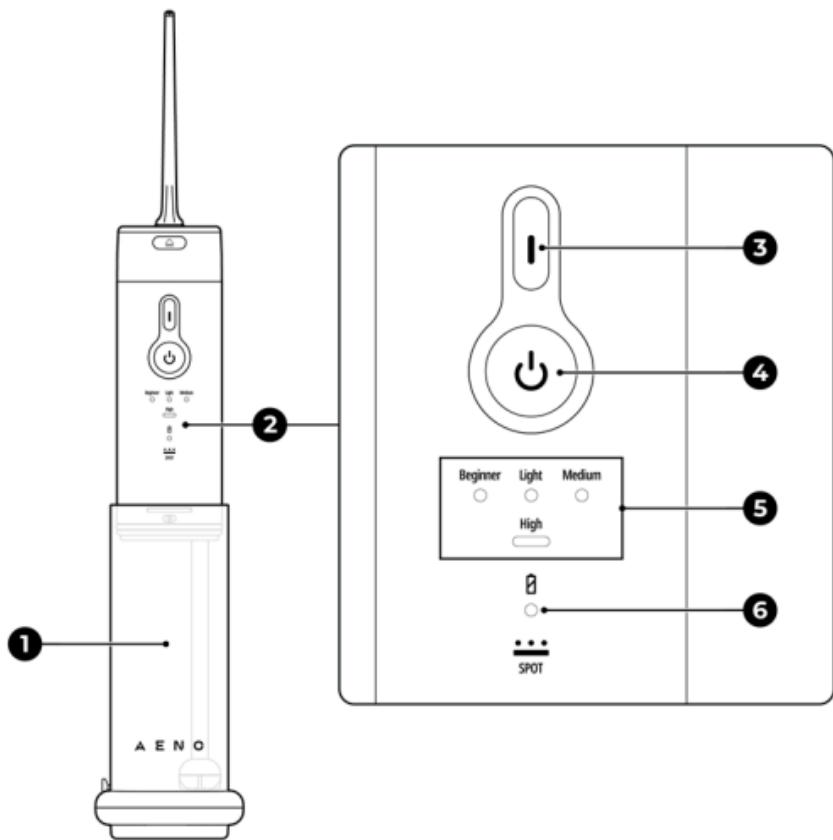
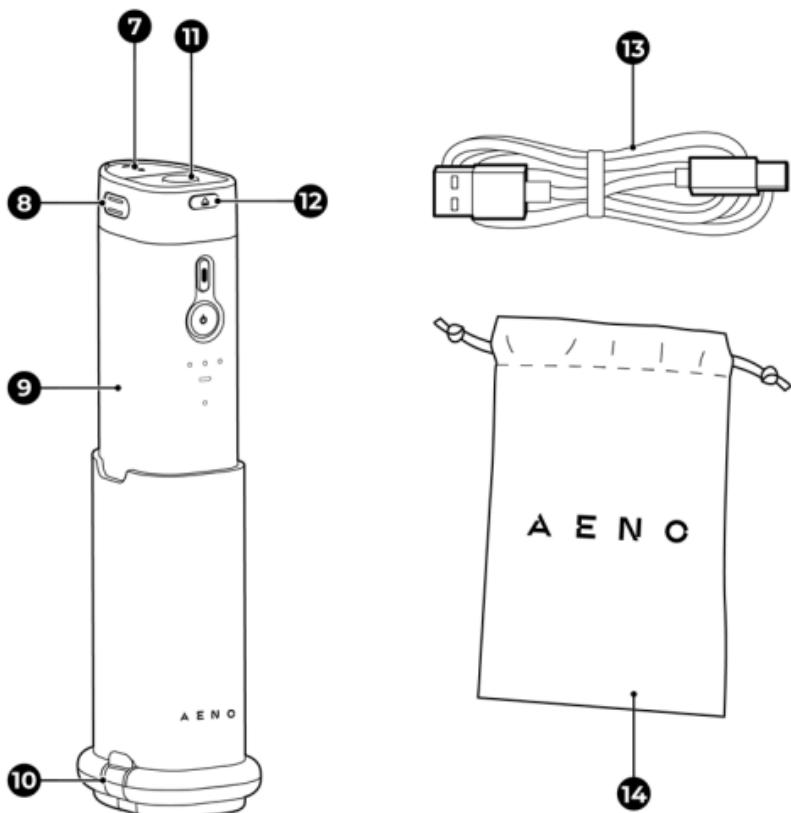
Dental Irrigator

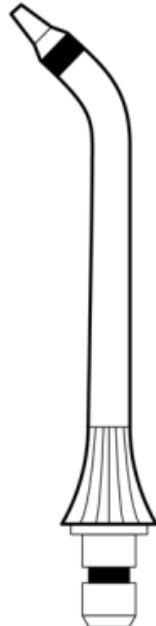
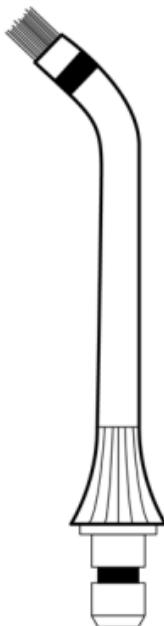
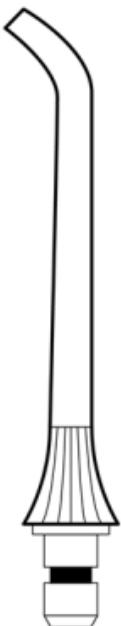
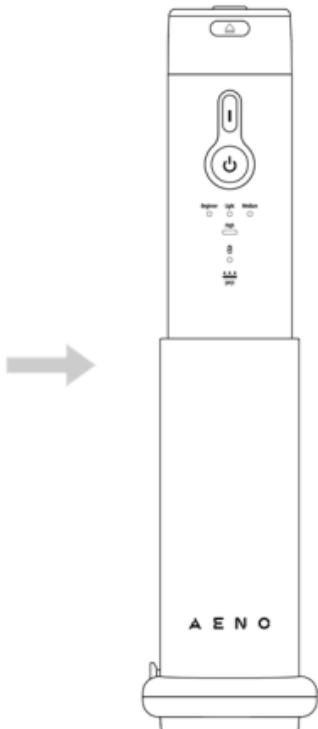
DI2

(Model: ADI0002)

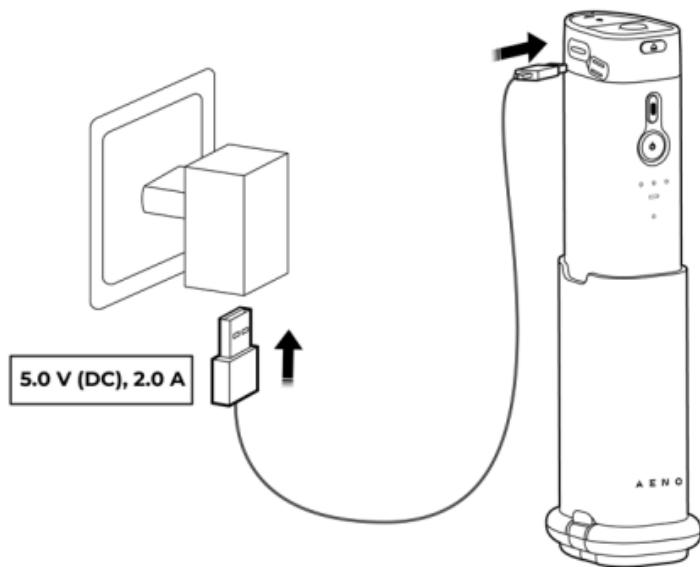
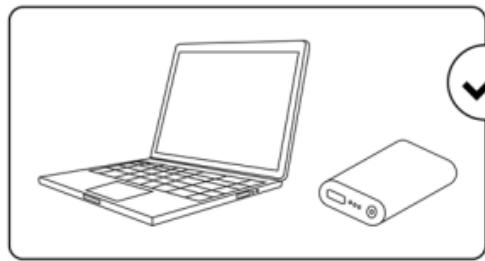
v 3.1.1 | 10.03.2025

ENG	6	KAZ	53
ARA	9	LAV	57
BOS	11	LIT	60
BUL	14	NLD	63
CES	18	POL	67
DEU	21	POR	71
ELL	25	RON	75
EST	28	RUS	78
FRA	31	SLK	82
HRV	35	SLV	85
HUN	38	SPA	88
HYE	42	SRP	92
ITA	45	UKR	95
KAT	49	UZB	99

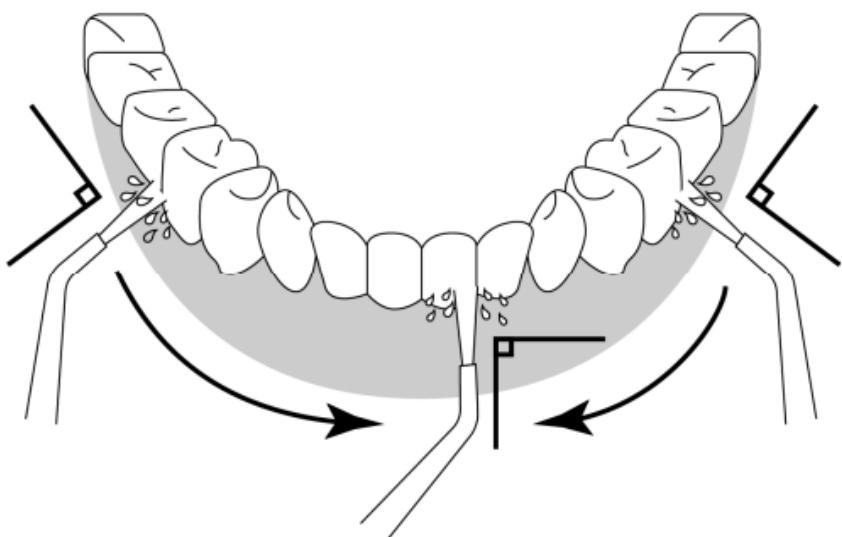
A**B**

C**15****16****17****D****E**

F

G

ENG The **AD10002** Dental Irrigator is designed to remove dental plaque and food debris from hard-to-reach areas of the mouth.

Technical Specifications

Power input: 5.0 V (DC), 2.0 A, 10.0 W (max.). Battery: non-removable, lithium-ion; 2,500 mAh, 3.7 V. Water pressure: 345–965 kPa. Pulsation: 1,900 pulses/min. Water reservoir: 200 ml. Full battery charge time: up to 3 h. Operating time on a single charge: at least 200 min. Operating time on a full reservoir: at least 40 s. Ingress protection rating: IPX7. Noise level (at a distance of 0.3 m): ≥ 68 dB. Control: manual. Color: White. USB (A) / USB (C) cable length: 0.6 m. Material: ABS plastic. Collapsed dimensions (L×W×H): 73×53×136 mm. Operating conditions: temp. +5...+40 °C, RH ≤ 75 %. Storage conditions: temp. 0...+45 °C, RH ≤ 75 %.

Package Contents (see Fig. **B, C**)

AD10002 Dental Irrigator, removable nozzles (3 pcs.), USB (A) / USB (C) cable, travel bag, Quick Start Guide.

Note. A power adapter is not included. The adapter output parameters must be 5.0 V (DC), 2.0 A.

Device Elements and Accessories (see Fig. **A, B, C**)

1 – water reservoir, 2 – control panel, 3 – operating mode selection button, 4 – power button, 5 – operating mode indicator, 6 – battery charge indicator, 7 – water reservoir fill door, 8 – USB (C) port with a cover, 9 – main body, 10 – water drain hole with a plug, 11 – opening for the nozzle, 12 – nozzle release button, 13 – USB (A) / USB (C) cable, 14 – travel bag, 15 – standard nozzle, 16 – orthodontic nozzle, 17 – periodontal nozzle.

Limitations and Warnings

This device is not intended for use by children under the age of 6 or by persons with reduced physical, mental or intellectual abilities unless they have sufficient experience and knowledge of how to use the device and are supervised by a person responsible for their safety. Do not let children play with the device and accessories. If you have oral cavity diseases, you should consult your dentist before using the device. Do not use a damaged device or accessories. If the device has been transported or stored at low temperatures, leave it at room temperature for 2 hours before use. Do not drop or throw the device. Clean the device only after turning it off and disconnecting it from the power supply. Do not use chemical or corrosive cleaning agents, abrasive pastes, agents that contain acids or solvents, or metal sponges to clean the device. Do not wash the device or its accessories in a dishwasher. Only use original accessories. Replace nozzles at least once every 3 months. Do not point the outlet of the nozzle into your eyes while the device is in use. If you are not using the device for a long time, charge it every 2–3 months to prevent deep discharge of the battery. For a detailed description of the device, its modes and features, see the full User Manual available at aeno.com/documents.

Intended Purpose of Nozzles

- Standard nozzle (15) is designed to remove plaque and food debris from between the teeth;
- Orthodontic nozzle (16) is designed to clean dental braces;
- Periodontal nozzle (17) is designed to clean and rinse periodontal pockets.

If you experience any issues or difficulties while using your AENO™ device, please email support at support@aeno.com or chat with us at aeno.com/service-and-warranty. Our specialists will help you solve the problem, and you will not have to spend time and effort visiting the store.

Button Functions

Name	Action	Function
Operating mode selection button (3)	Press once	Select operating mode
Power button (4)	Press once	Turn the device on/off
Nozzle release button (12)	Press and hold	Release the nozzle to remove it

Battery Charge Indication

Color	Status	Description
Green	On	The device is connected to a power supply and fully charged
		The device is turned on. The charge level is > 10 %
Red	Blinking	The device is connected to a power supply. The device is charging
		The device is turned on. The charge level is < 10 %
		The device's battery is fully discharged, the device does not turn on
-	Off	The device is turned off or not connected to a power supply

Operating Mode Indication

When the device is running in Beginner, Light, Medium or High modes, the indicator located on the control panel (2) under the name of the selected mode is solid white.

When the device is running in Spot mode, all 4 operating mode indicators on the control panel are solid white.

Operating Modes of the Device

- Beginner: low water pressure. Suitable for cleaning sensitive teeth and gums, as well as for first-time users.
- Light: moderate water pressure. Suitable for gentle cleaning. Recommended as the next mode after Beginner for a smooth transition to modes with a stronger water pressure.
- Medium: standard water pressure. Suitable for daily cleaning.
- High: high water pressure. Suitable for intensive cleaning. Recommended for experienced users.
- Spot: thorough cleaning of the oral cavity, removal of food debris from interdental spaces and prevention of tartar formation.

Note. In Spot mode, the device alternately ejects a jet of water for 2 seconds and then stops the water ejection for 2 seconds.

Preparing the Device for Use

1. Carefully open the package and take out the device and its accessories. Remove the device from the travel bag (14).
2. Pull the main body (9) all the way up (see Fig. D).
3. Detach the reservoir (1) as shown in Fig. E. Rinse it and the nozzles under running water, then wipe them with a soft, dry cloth.

Note. The reservoir locks in place tightly. You need to apply force to detach it.

4. Reattach the reservoir and put the device on a flat, stable surface.
5. Remove the cover from the USB (C) port (8) and connect the device to a power supply * using the USB (A) / USB (C) cable (13) (see Fig. F).
6. Fully charge the device (takes at least 3 hours). Once the battery is fully charged, the battery charge indicator (6) will turn green.
7. Disconnect the device from the power supply and insert the cover into the USB (C) port of the device.

* You can use mains electricity (a power adapter is required), desktop computers and laptops that have a USB (A) port, or portable chargers to charge the device.

Using the Device

1. Open the fill door (7) and fill the reservoir with water.

ATTENTION! Only use purified water without adding mouthwashes and other liquids or substances.

2. Close the fill door.

3. Choose the desired nozzle (see "Intended Purpose of Nozzles") and insert it into the opening (11) until it clicks into place.

Note. The nozzles can rotate 360°.

4. Set the device to the desired mode using the operating mode selection button (see "Operating Modes of the Device").

5. Point the outlet of the nozzle into the oral cavity and turn the device on by pressing the power button.

6. Carry out oral hygiene by moving the outlet of the nozzle along the edge of the gum as shown in Figure G.

Note. If the reservoir runs out of water during use, turn off the device by pressing the power button once and refill the reservoir as described in Steps 1 and 2 above. After that, continue performing your oral hygiene.

7. Turn off the device by pressing the power button once.

8. Remove the plug (10) and drain the remaining water.

9. Clean the device (see "Daily care").

Cleaning and Maintenance

Daily care. Direct the outlet of the nozzle from yourself and turn the device on by pressing the power button once. Wait for the device to stop ejecting water completely. After that, turn the device off by pressing the power button again. Press the nozzle release button and pull the nozzle up to remove it. Detach the reservoir as shown in Figure E. Rinse the nozzle and reservoir under running water, then dry them thoroughly. Reattach the reservoir and nozzle.

Total care. The reservoir is to be cleaned thoroughly once a week. Detach the reservoir as shown in Figure E. Wash it outside and inside under running water with a soft brush or sponge, without using cleaning agents, then dry it thoroughly. Reattach the reservoir.

Troubleshooting

The device does not turn on. Possible cause: the battery is discharged. Solution: charge the device's battery (charging time is 3 hours minimum).

The device does not eject water. Possible causes: the reservoir is out of water; the outlet of the nozzle is clogged. Solution: fill the reservoir with water; clean the outlet of the nozzle.

Insufficient water pressure. Possible cause: the wrong operating mode has been selected. Solution: set the required operation mode using the mode selection button.

ATTENTION! If none of the suggested solutions helped you eliminate the problem, contact the supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself.

ASBISC reserves the right to modify the device and to revise and change this document without prior notice to users.

The warranty period and service life are 2 years from the date of retail sale of the device.

Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Tel: +48732080077, info@asbis.com, asbis.com. Made in China. All trademarks used herein are the property of their respective owners. The date of manufacture is marked on the package.

Up-to-date information about and a detailed description of the device, as well as connection instructions, certificates, and details about the companies that deal with quality and warranty claims are available at aeno.com/documents. This device complies with applicable safety requirements.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

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ENG Recycling information. These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility. You should contact your local household waste disposal service for details.

For devices with a non-removable battery: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowable temperature, then, under abnormal conditions, the device will stop charging. If the battery temperature drops below the minimum allowable temperature during charging, then, under abnormal conditions, the current will be limited.

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