A E N O


User Manual

Toothbrush DB5 / DB6




(Models: ADB0005,
ADB0006)

v 2.1.0 | 11.04.2025


ENG.....	7
ARA.....	11
BOS.....	14
BUL.....	18
CES.....	22
DEU.....	26
ELL.....	31
EST.....	35
FRA.....	39
HRV.....	43
HUN.....	47
HYE.....	51
ITA.....	55
KAT.....	59
KAZ.....	64
LAV.....	68
LIT.....	72
NLD.....	76
POL.....	80
POR.....	84
RON.....	88
RUS.....	92
SLK.....	96
SLV.....	100
SPA.....	104
SRP.....	108
UKR.....	112
UZB.....	116

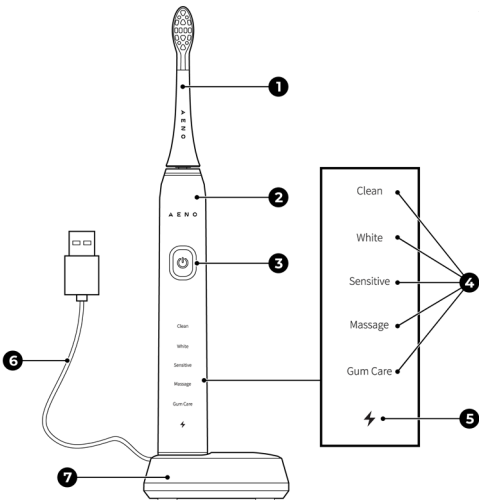
B

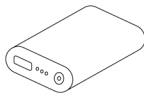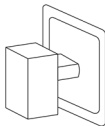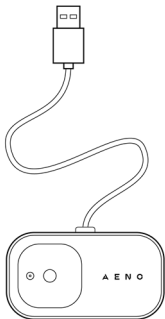C

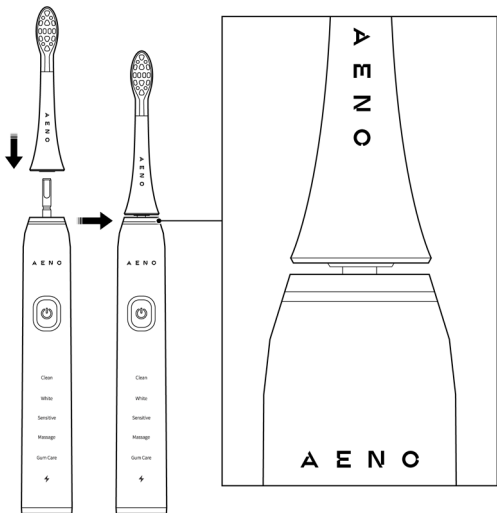D

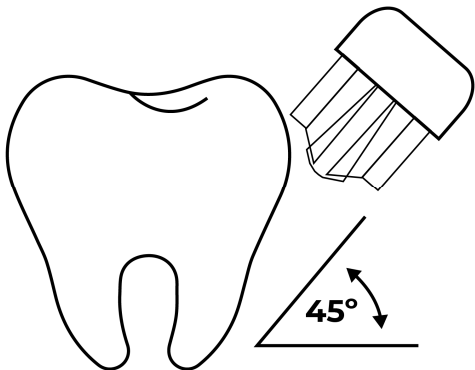ENG **AENO™** Toothbrush is intended for daily oral hygiene.

Models: ADB0005 (White), ADB0006 (Black).

Technical Specifications

Toothbrush. Input: 5.0 V (DC), 0.4 A, 2.0 W (max.). Battery: non-removable, lithium-ion, 750 mAh, 3.7 V, full charge time – up to 4.5 h. Operating time after a full charge: at least 130 min. Vibration intensity: up to 46,000 vibrations/minute. Number of operating modes: 5. Ingress protection rating: IPX7. Noise level (at a distance of 0.3 m): ≤ 65 dB. Control: manual. Dimensions of the toothbrush with the brush head (diameter × height): 27.6×246 mm. Casing material: ABS plastic. Bristle stiffness: medium. Service life of brush heads: 3 months. Weight (without brush heads): 104 g. Operating and storage conditions: temp. 0...+40 °C, RH 45–95 %.

Charging Station. Input: 5.0 V (DC), 1.0 A, 5.0 W (max.). Output: 5.0 V (DC), 0.6 A, 3.0 W (max.). Ingress protection rating: IPX7. Material: ABS plastic. USB (A) cable length: 1.2 m. Dimensions (L×W×H): 80×48×26 mm. Weight: 83.4 g.

Package Contents (see Fig. A)

AENO™ Toothbrush, replacement brush heads (2 pcs.), charging station with a USB (A) cable, Quick Start Guide.

Note. Power adapter is not included. The output parameters of the adapter must be 5.0 V (DC), 1.0 A.

Device Elements and Accessories (see Fig. A)

1 – replacement brush head, 2 – handle, 3 – button, 4 – operating mode indicators, 5 – battery charge indicator, 6 – USB (A) cable, 7 – charging station.

Limitations and Warnings

ATTENTION! If you have a pacemaker, oral diseases, or if you have just undergone oral surgery, consult your doctor before using the device.

If you experience any issues or difficulties while using your AENO™ device, please email support at **support@aeno.com** or chat with us at **aeno.com/service-and-warranty**. The specialists will help you solve the problem, and you will not have to waste time and effort visiting the store.

ATTENTION! To extend battery life, do not store the device on a charging station connected to a power source.

This device is not intended for use by children under the age of 6 or by persons with physical, mental or intellectual disabilities unless they have sufficient experience and knowledge of how to use the device and are supervised by a person responsible for their safety. Do not allow children to play with the device and its accessories. Do not drop or throw the device. Do not turn on or use a damaged device. If your gums begin to bleed after using the device and the bleeding continues for several days, contact your dentist. Use only the original charging station to charge the toothbrush. Make sure that the voltage and frequency specified in the technical documentation correspond to the power supply parameters. Each family member should use an individual brush head. Replace the brush head at least once every 3 months. Use only original replacement brush heads. If you are not using the device for a long time, charge it every 2-3 months to prevent deep discharge of the battery.

Statues of the Battery Charge Indicator

Color	Status	Description
Green	On	The device is placed on a charging station that is connected to a power source, the battery is fully charged
Red	Flickering *	The battery is charging
	Blinking **	Low battery level
–	Off	The device is turned on or placed on a charging station that is not connected to a power source
		The device is turned on, battery level is medium/high

* The indicator is on for 3 s, then off for 0.5 s.

** The indicator is on for 0.5 s, then off for 0.5 s.

Statues of the Operating Mode Indicators

The name of the current operating mode is highlighted on the device.

Button Functions

Action	Description
Press once (device is turned off)	Turn the device on

Action	Description
Press once (it has been less than 8 seconds since the device was turned on)	Switch operating modes
Press once (it has been more than 8 seconds since the device was turned on)	Turn off the device

Operating Modes

- “Clean” for daily brushing of the teeth;
- “White” for teeth whitening;
- “Sensitive” for brushing sensitive teeth;
- “Massage” for a thorough gum massage;
- “Gum Care” for a gentle gum massage.

Preparing the Device for Use

Charging. Carefully open the package and take out the device and its accessories. Clean them (see “Cleaning and Maintenance”). Place the charging station (7) on a dry, horizontal surface. Connect the station to a power supply (see Fig. B). Place the handle (2) on the charging station. The indicator (5) will start to flicker red. When the indicator turns green, the battery is fully charged and the device is ready for use. When charging is complete, disconnect the charging station from the power source.

Installing the Brush Head. Remove the handle from the charging station and install the brush head (1) as shown in Figure C.

Choosing the Operating Mode. Press the button (3) to turn on the device. Continue to press the button to switch the operating modes. You can change the mode within 8 seconds after turning the device on. Then press the button again to turn off the device.

Note. The last selected operating mode is stored in the device memory and will be activated the next time the device is turned on.

Operation

Wet the bristles of the brush head and apply the desired amount of a dental care product to the bristles. Point the bristles of the brush head toward the tooth surface at an angle of about 45° (see Fig. D). Turn on the device by

pressing the power button once. Brush your teeth without putting excessive pressure on the tooth enamel. Every 30 seconds, the device will notify you to change the brushing zone by briefly stopping the vibration. Change the brushing zone after each notification. After 2 minutes of operation, the device will turn off automatically.

Cleaning and Maintenance

ATTENTION! Turn off the device and disconnect the charging station from the power source before cleaning the device and accessories.

ATTENTION! Do not wash the device or its accessories in a dishwasher.

ATTENTION! Do not use chemical and aggressive detergents, abrasive pastes, products containing acids and solvents or steel wool scrubbers to clean the device and accessories.

Cleaning the Brush Head. Rinse the brush head, including the bristles, under running water.

Cleaning the Handle. Rinse the handle under running water or wipe with a soft cloth dampened with water. Then wipe dry.

Cleaning the charging station. Wipe with a soft cloth dampened with water. Then wipe dry.

Troubleshooting

The toothbrush does not turn on. Possible cause: the battery is discharged. Solution: charge the device (charging time – up to 4.5 hours).

The toothbrush is not charging. Possible causes: the charging station is not connected to a power source; the USB cable is damaged. Solution: connect the charging station to a power source (see Fig. B); check the integrity of the USB cable.

Weak vibration of the brush head. Possible cause: the place where the handle and the brush head connect is dirty. Solution: clean the handle and the brush head.

ATTENTION! If none of the suggested solutions helped you eliminate the problem, contact the supplier or an authorized service center. Do not disassemble the device or attempt to repair it yourself.

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ENG ASBISC reserves the right to modify the device and to revise and change this document without prior notice to users.

The warranty period and service life are 2 years from the date of retail sale of the device. Manufacturer: ASBISC Enterprises PLC, Iapetou, 1, Agios Athanasios, 4101, Limassol, Cyprus. Tel.: +48732080077, info@asbis.com, asbis.com. Made in China. All trademarks used herein are the property of their respective owners. The date of manufacture is marked on the package.

Up-to-date information about and a detailed description of the device, as well as connection instructions, certificates, and details about the companies that deal with quality and warranty claims are available at **aeno.com/documents**. This device complies with applicable safety requirements.

RoHS Directive. The device complies with the requirements of the Directive 2011/65/EU on the Restriction of Hazardous Substances (RoHS), including the requirements of the RoHS Directive 2015/863/EU.

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ENG Recycling information. These symbols indicate that you must follow the Waste Electrical and Electronic Equipment (WEEE) and Waste Battery and Battery Regulations when disposing of the device, its batteries and accumulators, and its electrical and electronic accessories. According to the rules, this equipment must be disposed of separately at the end of its service life. Do not dispose of the device, its batteries and accumulators, or its electrical and electronic accessories together with unsorted municipal waste, as this will harm the environment. To dispose of this equipment, it must be returned to the point of sale or turned in to a local recycling facility. You should contact your local household waste disposal service for details.

For devices with a non-removable battery: under normal or abnormal operating conditions, as well as in the event of a single malfunction, the charging voltage and current must not exceed their maximum permissible values. If the battery temperature exceeds the maximum allowable temperature, then, under abnormal conditions, the device will stop charging. If the battery temperature drops below the minimum allowable temperature during charging, then, under abnormal conditions, the current will be limited.

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